


bloomea®

PARIS



NeoCare

SKIN EXPERT

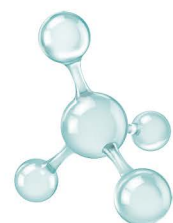

## WHAT IF YOUR EXPERTISE AND NEOCARE ENERGY WERE ONE ?

We have the answer: thanks to Tecar-Esthetique, **Bloomea** revolutionizes high-frequency current.

This latest-generation technology is the result of a French co-branding between **the leader** in high-frequency therapeutic currents, and **Bloomea**, skin treatment expert.

**NeoCare** meets the expectations of a very demanding male and female clientele by combining efficiency and comfort in the same treatment.

### YOUR EXPERTISE



# STIMULATE YOUR BODY NATURAL RESOURCES

**NeoCare** uses high-frequency currents between 300 and 1000 KHz, making it possible to work on 2 types of cells: **fibroblasts**, which are the cells of youth, and **adypocytes**, those of slimness.



OUR  
ENERGY

①

ACTS ON  
SKIN AGING

②

SCULPTS THE  
SILHOUETTE &  
ELIMINATE CELLULITE




## BIOLOGICAL ACTION REGENERATE

It is the action of reconstructing your tissues such as your dermis, for example. **BIONE0** will thus relaunch the intra and extracellular exchanges of the layers of your skin, which will allow it to be firmer and more elastic. Your fine wrinkles fade and your skin regains a natural radiance.


## THERMAL ACTION DYNAMISE

The complementary action of **THERMONE0** draws in depth the active ingredients (*nutrient, oxygen...*) which are directed on the treated area. The metabolic activity of the zone will allow a fast consumption of the active ingredients as well as their evacuation. But that's not all, the ability to restart blood and lymphatic circulation will help the treated area to destocking.

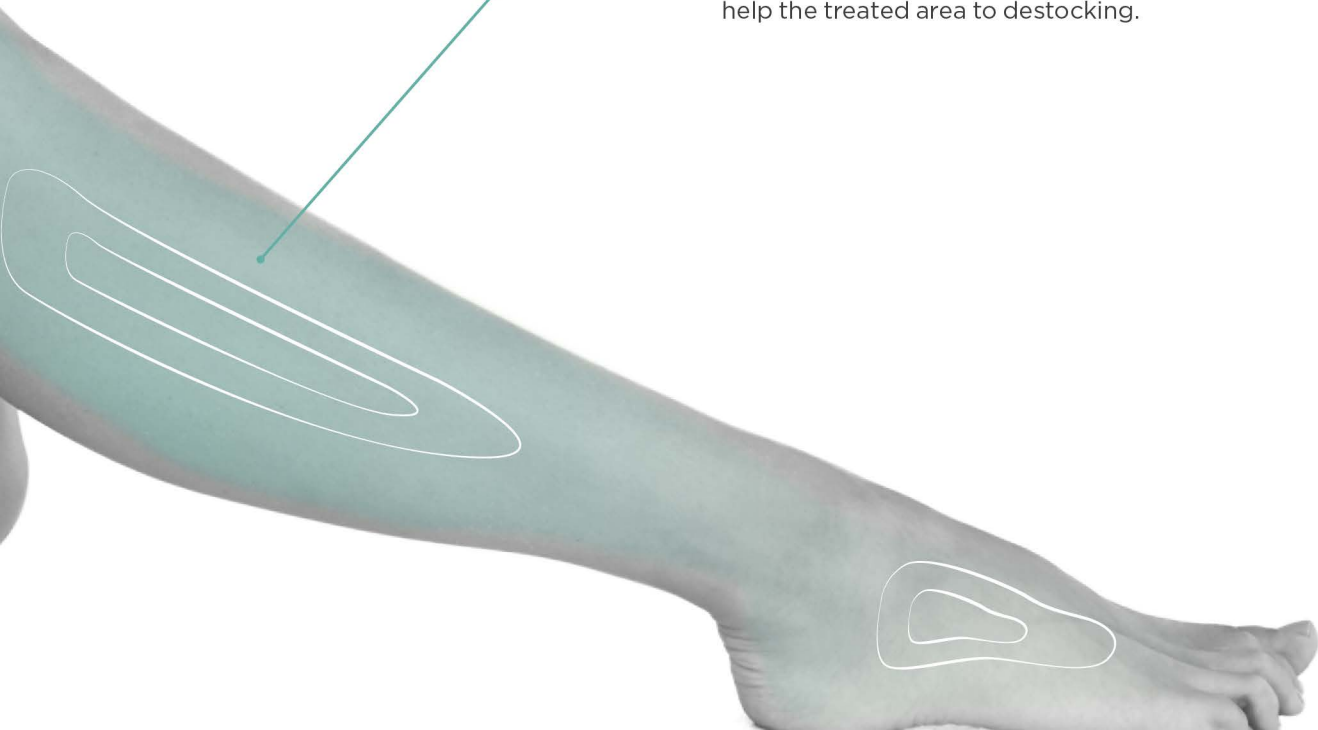



## NEO : SLIMMING BIO-BOOSTER

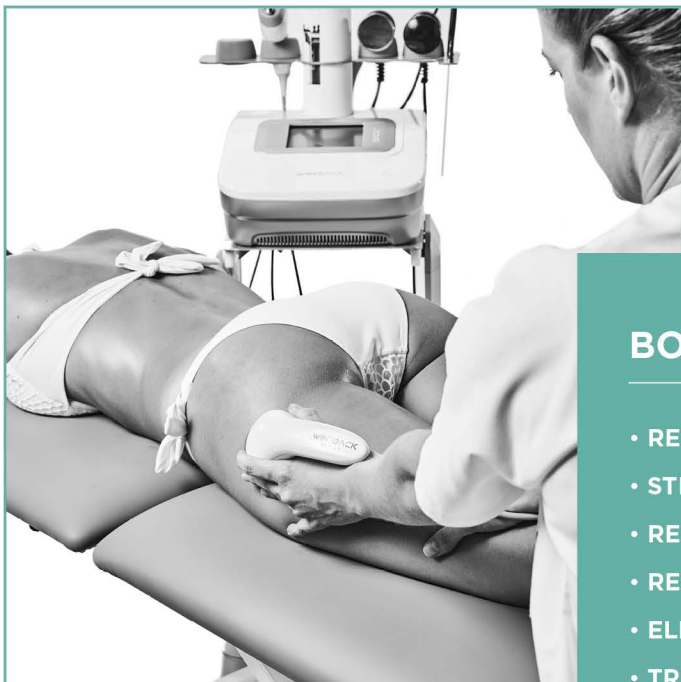

### BODY

---

- REDENSIFY THE SKIN
- STRENGTHEN ARMS
- REDUCE THE APPEARANCE OF CELLULITE
- RESHAPE CONTOURS
- ELIMINATE FLUID
- TREAT THE 3 TYPES OF CELLULITE






### **FIBROUS CELLULITE**

Present for a long time, very difficult to dislodge. Hard and painful to the touch.



### **ADIPOSE CELLULITE**

Often related to excess weight. Soft, thick and painless. Located on the belly, thighs and buttocks.

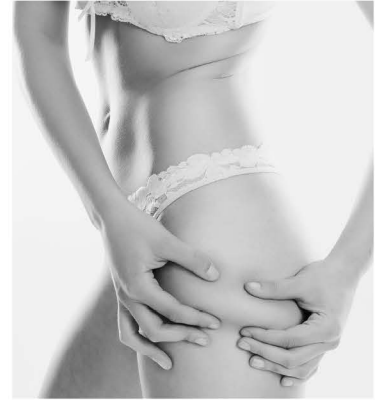

### **AQUEOUS CELLULITE**

Related to water retention due to poor venous and lymphatic circulation.



## NEO : YOUTH BIO-BOOSTER



### FACE & DECOLLETE

- **REFINE** THE SKIN TEXTURE
- **DRAIN** & RESTORE RADIANCE TO THE COMPLEXION
- **ELIMINATE** FLUID RETENTION (bags under eyes)
- **DIMINISH** WRINKLES AND IMPERFECTIONS (acne scars)
- **LIFT** FACE CONTOUR
- **REDUCE** DOUBLE CHIN
- **SMOOTH** THE NECKLINE
- **TREAT** THE 3 TYPES OF AGING



### SKELETONIZATION

Sunken cheeks / Deep wrinkles on the forehead / Marked eye contours / Sunken eye sockets / Creased lip contours



### IMPASTO

Double chin / Roundness of the face and cheeks / Weight gain



### SAGGING

Drooping eyelid / Sagging cheekbones / Marked nazogenian furrows / Jowls / Sagging oval of the face





## FREQUENCY PRINCIPLE

Scientific studies have shown that the frequency band between 300 and 1200 KHz induces an increase in intra- and extra-cellular exchanges and a diathermic effect on living tissue.


- **Monopolar current :**

The current passes between two electrodes. A mobile electrode (the handpiece) held by the practitioner and stationary electrode usually placed under the back.

- **Multipolar current :**

The current flows between the small electrodes on the handpiece.

Your **NeoCare** is equipped with 4 types of electrodes allowing different depths of action:

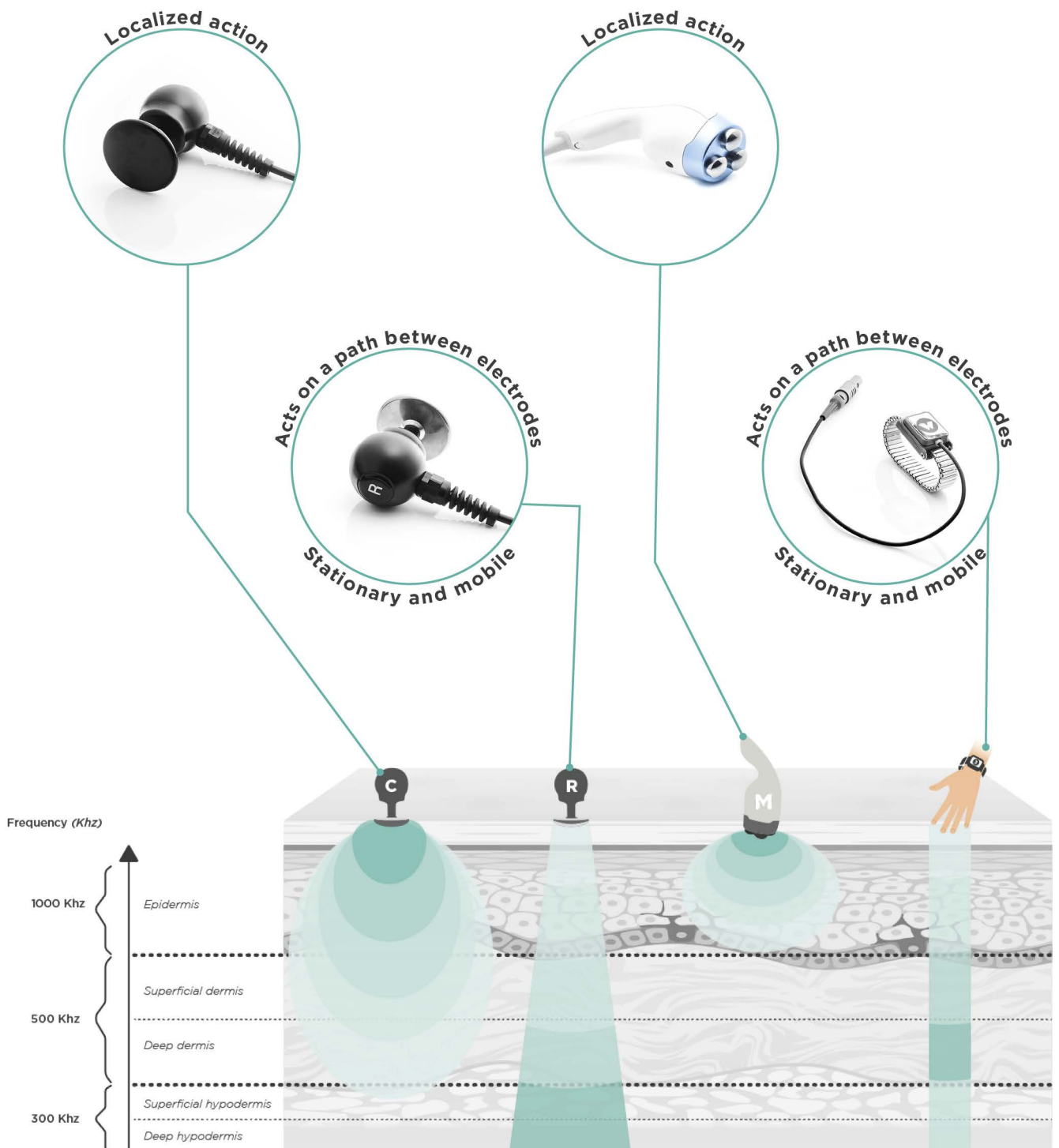

**NEOCARE INNOVATION:  
YOUR MANUAL KNOW-HOW  
IN THE SERVICE OF TECHNOLOGY.**

Thanks to a system of wristbands, NeoCare reconciles practical manual work and state of-the-art technology to deliver Neo energy to the deep dermis via a non-invasive current that is 100% natural for the human body.

The most delicate areas can be treated and the skin can be “shaped” at your convenience to treat skin sagging and the silhouette.



## USEFUL MODES

In order to fight against cellular addiction, NeoCare device uses 3 modes for an optimal and long-lasting result:



### DYNAMIC

Automatic frequency scanning between the superficial and deep layers of soft tissue. Simultaneous multilayer tension lifting & vascular recruitment.



CET



### LOW PULSE

Fine management on inflammatory areas / Acceleration of drainage techniques by playing on the lymph nodes / Optimization of the return to calm.



DEEP CET



### SUPER PULSE

Pulsed intensity peak allowing to work in hyperthermia while limiting the discomfort linked to the strong heat felt. Mainly used for fibrous, scar tissue to force areas with little vascularization to revascularize and promote collagen synthesis.



RET



WRISTBANDS

## FAST

The results are visible from the first session.

## LONG-LASTING

The desired effect is long-lasting by maintaining 1 session every month.

## COMFORTABLE

A session is completely noninvasive and painless.



BEAUTY

WELLNESS

# NeoCare

Additional applications

SCARS & STRETCHMARKS

DRAINAGE

DETOX





## NEOCARE PLUS

### NeoCare energy minus the weight

#### SPECIFICATIONS

---

Frequency	300 KHz, 500 KHz & 1 MHz
Modes	DYNAMIC - LOW PULSE - SUPER PULSE
Weight	19kg (device and trolley)
Dimensions	105x35x48 cm

#### ACCESSORIES

---

NeoCare 1.0	<b>MOBILE CET</b> : 3 electrodes of different sizes for focus or radial action <b>MOBILE RET</b> : 3 electrodes of different sizes for small and large areas <b>FIX</b> : Return plate and 10 FLEX adhesive plates <b>MyNeo</b> : Resistive wristband (RET)
NeoCare 4.0	<b>MIX3</b> (Body) - <b>MIX3</b> (Face)
NeoCare 5.0	2 adhesive plates RET, 2 RET wristbands Double resistive cable (RET)

*Equipment 2 years warranty - Accessories 1 year warranty  
Manufactured under standards ISO9001 - ISO14385 - KIWA 1984.*



*Time against  
the tide*

**Aesthetic Modeline  
Bloomea Paris**  
5 rue Poirier de Narçay  
75014 PARIS

**T. +33 (0) 9 86 34 45 73**  
contact@bloomea.fr  
[www.bloomea.fr](http://www.bloomea.fr)

**NeoCare**

SKIN EXPERT

